

# Varus / Valgus

Example 1



Example 2



Example 3



# With Your® Liner

Limb



3D model



Your® Liner



## Indication

The distal segment of the residual limb is mis-aligned with the load bearing axis of the leg. The below-knee segment of the stump is pointing inwards or outwards.



## Key features




Anatomically tailored

## Typical issues

### Off-the-shelf liner is over-stretched on external face

To conform with the stump, the standard liner will be stretched out on the external face, and can thus apply excessive pressure on this area. The stretched out cushioning layer might also offer insufficient cushioning, further exposing the underlying tissues.

The standard liner will also be distended on the internal face. This can lead to issues such as pinching of the skin, friction due to insufficient compression, or air pockets between liner and skin.

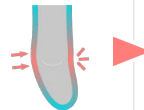

### Adapted shape and compressions

With its custom shape and adapted compressions, Your® Liner provides a close fit and adequate compressions all over the limb, eliminating the issue of excessive pressure, friction, air pockets and pinching of the skin.

The customisable thicknesses also allow to properly cushion sensitive areas, which can otherwise be exposed with absence of localised cushioning and overstretching of the external face



How can Your® Liner help ?